



Basic American Zenith® Z3 Specifications

Additional Features Include:

- High Quality UL 60601, IEC601-2-38 ,CAN/CSA C22.2 No. 601.1 compliant meeting the highest global safety standards for electrical medical equipment.
- Electrostatic powder coated antimicrobial sleep surface provides longer lasting, durable finish
- Grounded – Standard Hospital Grade
- Corner Mattress Retainers

U.S. Patents

- US 7,421,748 1
- Multiple Patents Pending

Technical Specifications:

- Mattress Deck Length..... 76"/80" optional 4" extender available
- Mattress Deck Width..... 35"/optional 42" extender
- Overall Bed Length..... 85.5" with bumper
- Overall Bed Width..... 39.9" with rails/assist devices
- Maximum back angle 70°
- Maximum knee angle 30°
- Maximum height 30"
- Minimum height 7.95"
- Horizontal translation less than 0.5"
- Safe Workload (including resident, support surface, head/foot boards, and accessories)500 lb

For more information on the Zenith Z3, custom options and available accessories, please visit us online at www.grahamfield.com or call Basic American today.



BASIC AMERICAN
MEDICAL PRODUCTS
800.554.9215

Specifications are subject to change without notice.
Basic American Medical Products is a registered trademark, SilverSolution and Zenith are trademarks of GF Health Products, Inc.
© 2009 GF Health Products, Inc. All rights reserved. GF0900152RevA09



The *Zenith Z3*, developed with both resident and caregiver in mind, is raising the bar in Long Term Care beds. While incorporating many of the same top quality features you've come to expect with the Zenith beds, the new *Zenith Z3* now offers unprecedented, advanced positioning speed that is faster than most of the


competition - thus saving time and reducing staff occupational injury when adjusting beds for resident care. With its superior design and quiet operation, the new *Zenith Z3* provides **speed, safety, and the solution** from a name you've come to trust.





Z3 Features:

- A) **Patent pending Dual Tie-Bar** design keeps bed level and balanced at all times providing a more comfortable, restful sleep
- B) **Unprecedented positioning speed** allows for bed positioning from a low 7.9" to 30" in as little as 26 seconds
- C) **Strong Supportive grid deck** is superior in frame design and offers a "no-sag" mattress support surface
- D) **Virtually 100%** vertical movement allows bed to take up less room space
- E) **Industry leading Z-lok foot caster base** with easy access one step red-green pedal locking system to maximize security, mobility, and when locked can be seen from a distance
- F) **360 Degree Mobility at ANY height** – can move from room to room at any height even the lowest.
- G) **Quiet, Smooth operation** by using 4 DC motors.
- H) **Adjustable Leg lift system** comes standard with Z3 for improving comfort for residents with edema and poor leg circulation
- I) **Wall Saver Bumper** comes standard and eliminates the need for rollers that damage walls
- J) **Customized headboards and footboards** are available in a variety of styles and finishes to compliment any room décor.

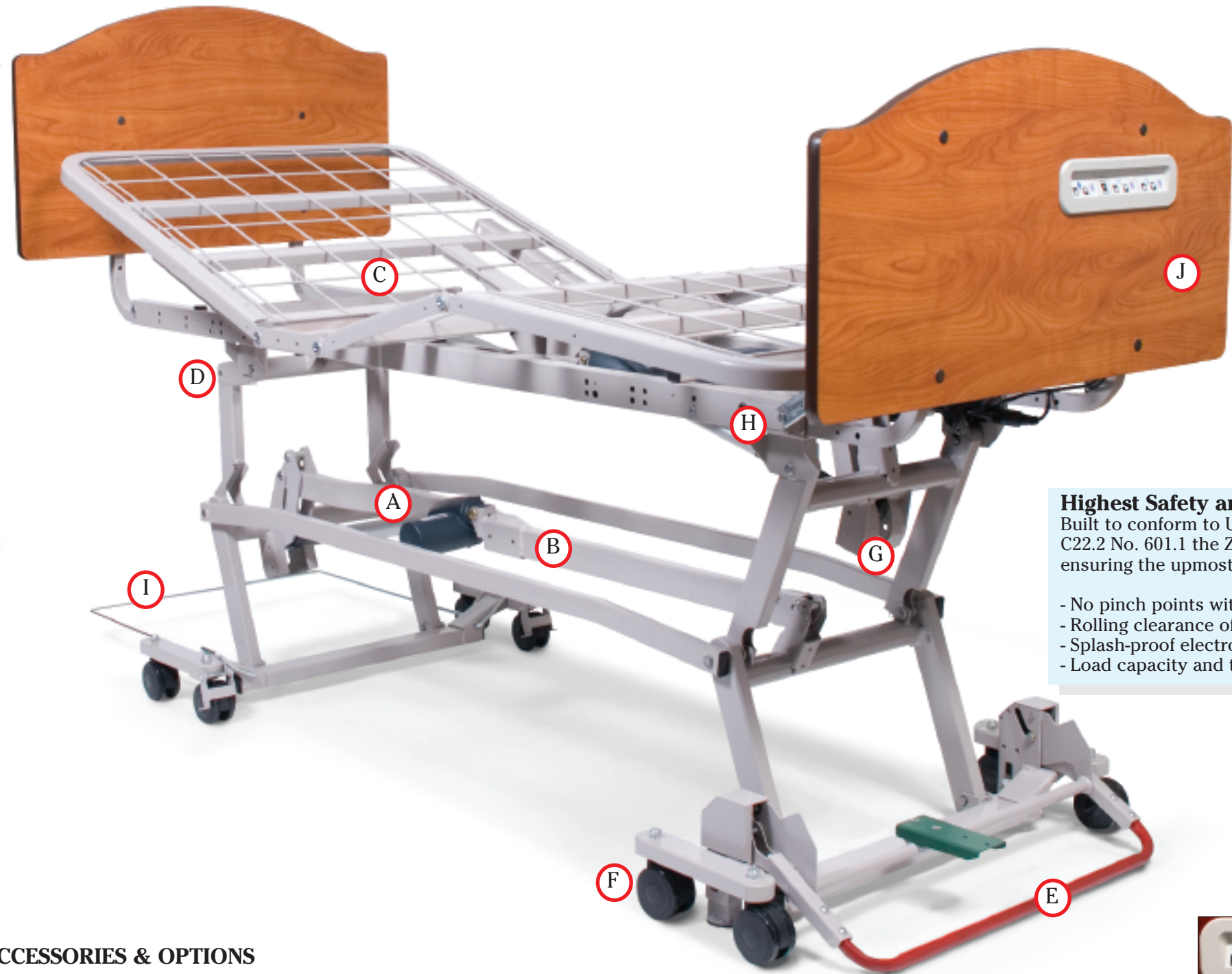

Also available with 42" Wide Deck option!

Highest Safety and Quality Standards

Built to conform to UL 60601-1, IEC601-2-38, and CAN/CSA C22.2 No. 601.1 the Zenith Z3 offers the following ensuring the upmost resident and staff safety:

- No pinch points within 200mm of bed perimeter
- Rolling clearance of 20mm over thresholds
- Splash-proof electronics
- Load capacity and tip-stability testing

And For Infection Control... Introducing SilverSolutions!

- Superior, high efficacy antimicrobial that protects all painted surfaces
- Antimicrobial kills more than 99.9% of odor and stain causing bacteria
- Lasts for the expected life of the bed
- Non-allergenic and safe for the environment
- High concentration of silver ions gives high efficacy
- Requires no special handling
- A cost effective solution that is built in to last

ACCESSORIES & OPTIONS

Extension Kit

Transform your bed to 84" improving inventory management of beds. Extension accommodates for a more broad range of resident heights.



Trapeze Adapter

Trapeze adapter allows for easy attachment of trapeze system to headboard with easy reach for residents.


Pendant Hand Control

Easy operation contoured Hand Control Pendant is shock-resistant, water-proof and durable. Holster conveniently secures to bedrails and assist devices for easy access.



ZFast

With *ZFast*, average speed for positioning is 1"/second at full 500 lb capacity. *ZFast* can reduce the chance for staff injury by reducing any bending and physical strain and improve resident care.

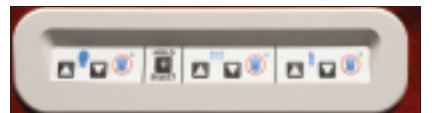

Staff Control

Optional Embedded Staff Control allows staff to operate bed functions from the footboard. Staff can restrict access to the resident by locking out control from the patient pendant. Two Stage locking helps prevent unwanted adjustment.

Half Length Assist Device and Assist Bar

Half Length Assist Device and assist bar provide a sturdy and secure handhold to assist residents in and out of bed. Assist devices are removable without tools and meet FDA Entrapment Guidelines.

