

Clinical Care Recliner

- Welded steel frame is constructed of round tubing resulting in improved structural integrity of the frame.
- Improved wider and softer armrest for patient comfort.
- Softer seat and back package for improved patient comfort. Varied density seat back and redesigned seating position provides excellent comfort and support during long term sitting.
- Easy to activate gas cylinder mechanism for quick Trendelenburg activation and improved lift.
- Durable, high impact side tables sit flush with chair arm to conveniently hold prep materials. Easy to actuate by caregiver.
- Easily removed “pop-out” side panels enable easy access to recliner for cleaning and maintenance.
- Three-position, plus Trendelenburg.
- Four Tente® swivel casters facilitate transport and maneuverability in confined areas. Left and right rear top-locking casters provide positive lock. Directional wheel lock on right rear caster reduces “crabbing” of recliner and enables the chair to be moved down hallways with less effort.
- Weight Bearing Capacity: 290 Pounds.

Specs

Overall Height.....	47"	Seat Depth.....	21"
Overall Width.....	31"	Depth Fully Reclined	
Overall Depth.....	35"	In Trendelenburg	74"
Width Between Arms	22"	Width With Both Side Tables	
Arm Height From Floor.....	25.5"	in Upright Position	51"
Arm Height From Seat.....	8"	Shipping Weight	127 lb
Seat Height.....	22"	Max. Weight Capacity	290 lb
Seat Width.....	21.5"		



PERMABLOK³
ADVANCED VINYL PROTECTION



Featuring specially-contoured high-density foam, the new back profile enhances patient comfort and support.



Easily removed “pop-out” side panels enable easy access to recliner for cleaning and maintenance.

Features **LIQUICELL[®]** Technology

LiquiCell is an ultra-thin, liquid-filled interface specifically designed and proven to significantly reduce shear stress, aiding in the prevention of pressure sores and general discomfort.



Infinite Position Back Design


Item #	Item # with Heat & Massage	Color
Patriot Plus Fabric		
FR577RG427	FR577RGH427	Blue Ridge
FR577RG432	FR577RGH432	Imperial Blue
FR577RG454	FR577RGH454	Royal Blue
FR577RG851	FR577RGH851	Doe Skin
FR577RG857	FR577RGH857	Jade
FR577RG863	FR577RGH863	Rosewood
FR577RG8504	FR577RGH8504	Sand
FR577RG8567	FR577RGH8567	Dove
FR577RG8602	FR577RGH8602	Oatmeal
FR577RG8604	FR577RGH8604	Concrete
FR577RG8561	FR577RGH8561	Mulberry
FR577RG8611	FR577RGH8611	Terra Cotta
FR577RG8615	FR577RGH8615	Berry
FR577RG8533	FR577RGH8533	Port
FR577RG8622	FR577RGH8622	Gray Moss
FR577RG8624	FR577RGH8624	Granny Smith
FR577RG8626	FR577RGH8626	Eucalyptus
FR577RG8584	FR577RGH8584	Taupe
FR577RG8516	FR577RGH8516	Chestnut
FR577RG8628	FR577RGH8628	Sea

Savoy Fabric		
FR577RG401	FR577RGH401	Vintage Gold
FR577RG409	FR577RGH409	Warm Taupe
FR577RG9711	FR577RGH9711	Tallow
FR577RG9754	FR577RGH9754	Sisal
FR577RG9757	FR577RGH9757	Tumeric
FR577RG9701	FR577RGH9701	Bay Leaf
FR577RG9722	FR577RGH9722	Desert Sage
FR577RG9715	FR577RGH9715	Eucalyptus
FR577RG9764	FR577RGH9764	Steel Grey
FR577RG9728	FR577RGH9728	Navy

Renaissance Fabric		
FR577RG9200	FR577RGH9200	Chalk
FR577RG9205	FR577RGH9205	Gypsum
FR577RG9203	FR577RGH9203	Pyrite
FR577RG9202	FR577RGH9202	Periwinkle
FR577RG9201	FR577RGH9201	Cobblestone
FR577RG9210	FR577RGH9210	Cypress
FR577RG9207	FR577RGH9207	Ice Blue
FR577RG9214	FR577RGH9214	Steel Blue

Silvertex with SilverGUARD Fabric		
FR577RG8808	FR577RGH8808	Champagne
FR577RG8809	FR577RGH8809	Taupe
FR577RG8810	FR577RGH8810	Chestnut
FR577RG8814	FR577RGH8814	Umber
FR577RG8815	FR577RGH8815	Borscht
FR577RG8818	FR577RGH8818	Sage
FR577RG8819	FR577RGH8819	Celery
FR577RG8801	FR577RGH8801	Marine Blue
FR577RG8822	FR577RGH8822	Storm
FR577RG8805	FR577RGH8805	Jet

Refer to GF Health Products, Inc. Color Card

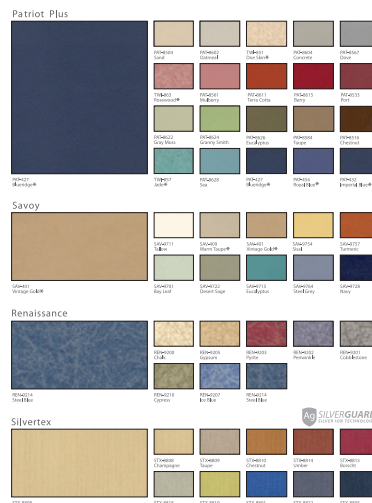

FR Prefix designates that the seating meets CA133 Fire Code Specifications.

PERMABLOK³® PROTECTION!


ADVANCED VINYL PROTECTION

All Lumex Specialty Healthcare Seating Upholstery features PERMABLOK3 Advanced Vinyl Protective Coating. Engineered to create a tough effective barrier against the three biggest problems encountered in healthcare environments – germs, abrasion, and stains. You can't see it or feel it, yet the protection it offers results in extended product performance and lasting beauty with a minimum of care. It guards against the surface growth of fungus, mold and mildew spores, which can cause unpleasant odors, unsightly pink and black stains, and even allergic reactions. All Lumex Upholstery Selections not only keep their good looks but also remain highly resistant to deterioration.

HOW Ag SILVERGUARD[®] SILVER ION TECHNOLOGY WORKS!

Silver Ions destroy microorganisms by attaching to the bacteria's cell wall to disrupt protein metabolism within the cell. These are slowly released from their inorganic matrix in the presence of moisture. Simply stated—wipe down the upholstery with a damp cloth to activate the Silver Ions!

Silver encapsulates MRSA bacterium and it dies. SILVERGUARD is released as needed to effectively control bacterial growth, while leaving sufficient quantities of SILVERGUARD to retain antimicrobial protection.

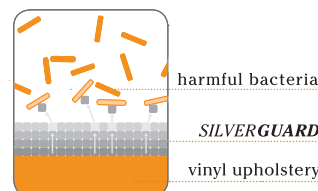

GRAHAM-FIELD

www.grahamfield.com

© 2010 GF Health Products, Inc.
 Atlanta, Georgia 30360 GF0900050RevB10