

The Liberty Bed provides you with the key features you want at a very affordable price point.

Strong, no-sag grid deck sleep surface with sturdy cross braces. Safe work load: 450 pounds*

• Height travel range: 9" to 24.4".With a standard 6" mattress, this positions the bed at an ideal height for the caregiver and for transfers.



• One bed can be used for 76" and 80" lengths using the adjustable brackets to meet space limitations and improve inventory management. No extension kit is necessary.





- Resident friendly, quiet, smooth operation keeps bed level and balanced at all times using three DC motors.
- 360 degree mobility at <u>ANY</u> height can move from room to room at any height, and move over thresholds at even the lowest position.
- Mattress retainer is standard.
- **Reinforced**, **robust base** is designed for support yet provides sufficient underbed clearance, ideal for using overbed

tables and patient lifts.

 Easy assembly requires only a screwdriver.



Keyed Wall Saver is standard – easily moves over thresholds during transport, yet automatically falls back into place to protect your walls. Wallsaver adjusts easily when modifiying bed from 76" to 80" and when using a trapeze.

Pendant Hand Control

Easy operation contoured Hand Control Pendant is shock-resistant, water-proof and durable. Holster



conveniently secures to bedrails and assist devices for easy access.



ACCESSORIES & OPTIONS



Optional Trapeze Adapter and **Trapeze** attaches easily to headboard to position trapeze within easy reach of residents.

ZA83000 (Adapter for use with ZA781000 Trapeze Bar)

Optional Liberty Lock One-step red/green pedal locking system maximizes security. Locked position can be seen from a distance. The Liberty Brake can be



With Pedal Lock Four Locking Casters secure the



added as a retrofit.

bed at any height. LBT7174 With 4 Locking Casters



Highest Safety and Quality Standards

Built to conform to UL 60601-1 and IEC601-2-38, CAN/CSA C22.2 No. 601.1 the Liberty Bed offers the following ensuring the utmost resident and staff safety:

- No pinch points within 200mm of bed perimeter.
- Rolling clearance of 20mm over thresholds.
- Splash-proof electronics.
- Load capacity and tip-stability testing.

And For Infection Control... Introducing SilverSolutions! • Superior, high efficacy that protects painted surfaces

- Antimicrobial kills more than 99.9% of odor and stain causing bacteria.
- Lasts for the expected life of the bed.
- Non-allergenic and safe for the environment.
- High concentration of silver ions gives high efficacy.
- Requires no special handling.
- A cost-effective solution that is built-in to last.



Optional Extension kit

Extends sleeping surface from 80" to 84". LBT83300

Optional Half Length Assist Device and Assist Bar

Half length assist device and assist bar provide a sturdy and secure handhold to assist residents in and out of bed. Assist devices are removable without tools and meet FDA Entrapment Guidelines.






LBT80300



ZA82800 (head end, pair) ŽA82900 (foot end, pair)



ZA85500

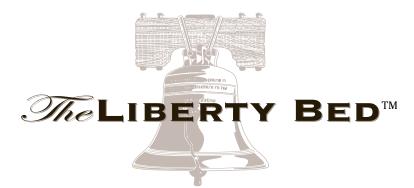

Basic American Liberty Bed[™] Specifications

Additional Features Include: Electrostatic powder-coated antimicrobial bed surface provides longer lasting, durable finish

LBT7174 The Liberty Bed with four locking castersLBTL7174 The Liberty Bed with pedal lock

Specifications:

Height Travel Range	9.9" to 24.4"
Safe Work Load*	
Mattress Deck Length	. 76"or 80"
Mattress Deck Width	35"
Overall Bed Length	86" with bumper
Overall Bed Width	39.75" with rails/assist
	devices
Maximum back angle	. 70°
Maximum knee angle	. 30°

Proudly Made in America



*including resident, support surface, head/foot boards, and accessories.

For more information on the Liberty Bed, custom options and available accessories, please visit us online at **www.grahamfield.com** or call Basic American today.



800.554.9215 Fax 800.524.4113 View these and more products at www.grahamfield.com

Basic American Medical Products is not responsible for typographical errors. All illustrations, specifications, packaging and warranties contained in this literature are based on the latest product information available at the time of printing. The most current product information can be found online at www.grahamfield.com. SilverSolutions and The Liberty Bed are trademarks of GF Health Products, Inc. GF1000040RevB12