

# HEALTH CARE SEATING



## HIGH PERFORMANCE FABRIC COLLECTIONS

**eco<sup>sense</sup> BY SPRADLING GREEN COLLECTION**

**GUARDIAN** eco<sup>sense</sup> BY SPRADLING

GUA-1325   <b>ARMOR</b>	GUA-1305   <b>CONVOY</b>	GUA-1323   <b>MOAT</b>	GUA-1319   <b>SIREN</b>
GUA-1324   <b>SENTRY</b>	GUA-1311   <b>FENCE</b>	GUA-1318   <b>STEEL</b>	GUA-1304   <b>LATCH</b>

**SAVOY** eco<sup>sense</sup> BY SPRADLING

SAV-9517   <b>EUCALYPTUS</b>	SAV-9524   <b>BLACK</b>	SAV-409   <b>WARM TAUPE</b>	SAV-9516   <b>INDIGO</b>

**CASSARA** eco<sup>sense</sup> BY SPRADLING

CAS-9611   <b>AMETHYST</b>	CAS-9613   <b>CORK</b>	CAS-9607   <b>PATINA</b>	CAS-9604   <b>INDIGO</b>

## PERFORMANCE COLLECTION

**DOLCE<sup>+</sup>**

DOL-6727   <b>METEOR</b>	DOL-6717   <b>MOSS</b>	DOL-6724   <b>JET</b>	DOL-6725   <b>INDIGO</b>

Some products may contain chemicals, identified by the State of California as a cause of cancer, birth defects or other reproductive harm at certain exposure levels. Under the State of California law known as "Prop 65", certain warnings to consumers in California are required where such products are sold. The presence of such chemicals in coated fabric products, even at extremely low concentrations, can lead to disputes over Prop 65 warning obligations that are very difficult to resolve.

If your company intends to introduce Spradling vinyl products into commerce in California, Spradling recommends that you take steps to confirm whether the warning requirements of California's "Prop 65" are applicable to the products in question and whether your company is complying with those warning requirements for the chemicals contained within these products.



Composition	NEO					
	Dolce	Patriot Plus	Renaissance	Guardian	Savoy	Cassara
PVC	■	■	■	■	■	■
PU	■					
TPE						■
100% Polyester Hi-Loft™ Backing			■		■	
Polyester / Cotton Backing		■				
Ultra Hi-Loft™ Polyester Backing				■		
100% Rayon	■					
100% Polyester						■
Weight (oz)	15	34	32	30	30	22
Roll Size (yd)	30	30	30	30	30	30
Roll Width (inch)	54	54	54	54	54	54

Codes*	NEO					
	DOL	PAT	REN	GUA	SAV	CAS
NFPA 260- Cover Fabric - Class 1	■	■	■	■	■	■
UFAC Fabric - Class 1	■	■	■	■	■	■
CAL TB 117-2013	■	■	■	■	■	■
FMVSS 302	■	■				■

Finishes and Treatments	NEO					
	DOL	PAT	REN	GUA	SAV	CAS
<b>PBG<sup>+</sup> PERMABLOK<sup>®</sup> PERMAGUARD<sup>®</sup></b> ADVANCED VINYL PROTECTION				■	■	
<b>SUREGUARD™</b>						■
Antibacterial Protection (1,000 hrs)				■		
AATCC TM 16-3 / NTC 1479 / CFFA-2a - Blue Wool (1,000 hrs)	■	■	■	■	■	■
QUV Stabilized Pigments (650 hrs)	■					
Abrasion - CFFA-1ASTM D-4157-02 (k)	100	100	100	100	100	100
Wyzenbeek Cotton Duck	#10	#10	#10	#10	#10	#10
Cold Crack (F)	-10°	-20°	-20°	-10°	-10°	-10°

Applications	NEO					
	DOL	PAT	REN	GUA	SAV	CAS
Corporate Interiors   Restaurant   Hotel	■	■	■	■	■	■
Casino   Sunroom   Residential	■	■	■	■	■	■
Healthcare	■	■	■	■	■	■
Marine Interior		■	■	■		
Exterior Approved				■		

published November 2023

**Color Variations** may occur due to monitor or printer settings. Swatches are digital reproductions. Perceived texture and coloring may vary on finished product.

**Fabric Availability** current at time of publication and subject to change without notice. Contact Customer Service for most current availability

\* This term and any corresponding data refer to the typical performance in the specific tests indicated and should not be construed to imply the behavior of this or any other material under actual fire conditions



## PREMIERE COLLECTION

### PATRIOT PLUS

NEO



PAT-432 | IMPERIAL BLUE



PAT-8602 | OATMEAL



PAT-8604 | CONCRETE



PAT-8567 | DOVE



PAT-454 | ROYAL BLUE




PAT-427 | BLUERIDGE



PAT-8584 | TAUPE



PAT-8516 | CHESTNUT

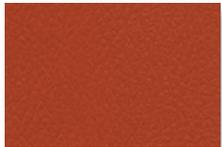

PAT-8611 | TERRA COTTA



PAT-8504 | SAND




PAT-8533 | PORT



PAT-8626 | EUCALYPTUS

### RENAISSANCE

NEO

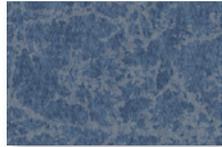


REN-9214 | STEEL BLUE




REN-9201 | COBBLESTONE




REN-9205 | GYPSUM

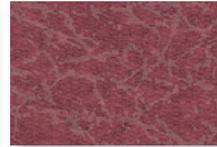


REN-9208 | WINEBERRY

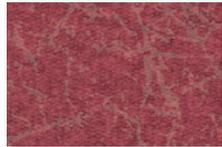

REN-9203 | PYRITE

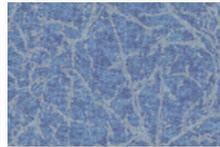

REN-9207 | ICE BLUE

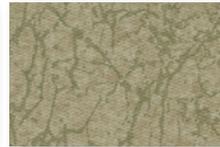

REN-9210 | CYPRESS



REN-9212 | BLUEJADE

## Care and Cleaning Guide†



We recommend cleaning at least once a week to extend the lifespan of the upholstery, which will prevent dirt and contamination buildup and permanent staining. In the event of a spill, wipe off any excess right away with a dry white cloth or fabric, but careful not to spread to remaining surface.



Clean the whole surface in a circular motion using liquid hand soap (pH neutral) and water (1 part soap to 9 parts water). Check the Approved Cleaners guide for recommended cleaning products. Wipe with a clean damp cloth to get rid of any excess soap. Repeat as needed and then let dry.



If there is any buildup that cannot be removed with soapy water, a mixture of 1 part isopropyl alcohol to 9 parts water can be used as a last resort. (Remember, using alcohol on a regular basis can affect the properties of the product).



Once a month, you can use a soft-bristle brush for a deeper cleaning. Remember, light colors need more care, more often.

### Disinfectants

In addition to cleaning, healthcare facilities maintain disinfecting regimens to reduce the spread of infections. There is a wide variety of available biocides that include bleaches, peroxides and quaternary ammonium products, among others. They should always be used at the recommended dilution, never in concentrated form. To prolong the life of coated fabrics, the disinfectants should be rinsed off after each use. It should be noted that misuse of disinfectants is the major source of surface deterioration. Spradling International, Inc.<sup>®</sup> continually evaluates cleaners/disinfectants, and a list of those tested can be found on our website.

Abrasion test results exceeding ACT Performance Guidelines are not an indicator of product lifespan. Multiple factors affect fabric durability and appearance retention.

### Note

The cleaning products recommended were used in a laboratory setting. Spradling Group<sup>®</sup> has no relationship with the manufacturers of those cleaning products.

Please visit [www.spradling.group](http://www.spradling.group) for the Polyurethane Fabric Cleaning Guide.



GF Health Products, Inc. is an  
intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.  
The most current and complete product information can be found on our website.

[www.grahamfield.com](http://www.grahamfield.com)



© 2015, GF Health Products, Inc. All Rights Reserved. ©2021 Spradling International, Inc.<sup>®</sup>  
PBG+<sup>®</sup>, PERMAGUARD<sup>®</sup> and PERMABLOK<sup>®</sup> are registered trademarks of Spradling Holdings, Inc.  
Lumex is a registered trademark of GF Health Products, Inc.

### For Color Transfer

Certain clothing and accessory dyes (such as those used on denim jeans) may migrate to lighter colors. This phenomenon is increased by humidity and temperature and is irreversible. Please check compatibility when using this product in combination with painted or varnished surfaces  
Color Variations

### Color Variations

may occur due to monitor or printer settings. Swatches are digital reproductions. Perceived texture and coloring may vary on finished product.

### Fabric Availability

current at time of publication and subject to change without notice. Contact Customer Service for most current availability.

† Care and Cleaning Guide excludes Dolce.