

★★★  
**PATRIOT LX**  
*Semi-Electric Homecare Bed*




- Robust steel grid deck sleep surface provides superior mattress support and easy to clean surface compared to traditional spring decks
- Head and foot sections can be raised or lowered simultaneously or independently with one hand control
- Recessed main frame offers maximum pinch point protection between grid deck and frame – reducing the risk of injury
- Lightweight design weighs only 118 lbs – lightest bed on the market
- Bed setup is so simple that it can save on setup time and requires minimal maintenance keeping your technicians time efficient
- Patent pending integrated HI-LO mechanism makes bed height adjustment easy and reduces potential loss of parts such as HI/LO rod when used as rentals
- Self contained crank for semi-electric that can be used at either end makes adjustment from low 10" to 20" easy and even allows for trendelenburg and slant positioning
- Patent pending rail slots in frame ensure quick and proper rail placement
- Single motor unit is designed for easy installation, requires less maintenance and provides the ultimate in quiet, smooth operation
- Attractive, “home-like” thermo fused bed ends offer a less institutional look and eliminate a bed component that often requires frequent repair and maintenance
- Features SilverSolutions- an anti-microbial technology built into the paint that protects against bacteria and odor
- During a power failure, the motor can lower both the head and foot section with just a single 9 volt battery
- 4" extension kit can be installed on head, foot or both to extend bed length to 84" or 88"

**Grid Deck offers the following advantages:**

- Much stronger frame
- Better mattress support
- No spring sag or bottoming out with heavier patients
- Easier to clean
- No connecting links or springs to lose

★★★  
**PATRIOT LX**  
*Semi-Electric Homecare Bed*

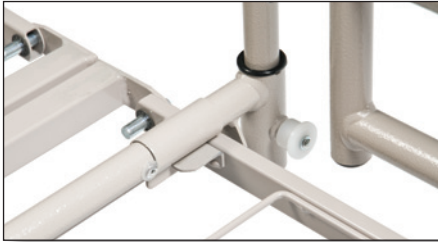

Rail Slot



Casters (optional)



Footboard

**Technical Specifications**

Total bed weight assembled: ..... 118 lb  
 Maximum patient weight..... 350 lbs EVENLY DISTRIBUTED  
 450 lb distributed total weight capacity  
 including patient, mattress, rails and accessories  
 Overall Bed Size:..... 84" L x 35" W  
 Lift Range:..... Head – 60, Knee – 30  
 Bed Height:..... Min. 10", Max. 20" (at deck level)

**Warranty:**

Bed Frame..... 5 Years  
 Electrical Components..... 2 Years  
 All Other Durable Components..... 1 Year

<b>US5000</b>	Semi-Electric Bed Only
<b>US5000-IFPKG</b>	Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Full Rails
<b>US5000-IFPKGHR</b>	Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Half Rails
<b>US5000-IFPKGQR</b>	Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Quarter Rails
<b>US5000-PKG</b>	Semi-Electric Bed with 1633-Innerspring Mattress and FDA Full Rails
<b>US5000-PKGHR</b>	Semi-Electric Bed with 1633-Innerspring Mattress and FDA Half Rails
<b>US5000-PKGQR</b>	Semi-Electric Bed with 1633-Innerspring Mattress and FDA Quarter Rails
<b>US5000-PKGC</b>	Semi-Electric Bed with 1633-Innerspring Mattress and Chrome Full Rails
<b>US5000-PKGNR</b>	Semi-Electric Bed with 1633-Innerspring Mattress and No Rails
<b>US5000-RFPKG</b>	Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Full Rails
<b>US5000-RFPKGHR</b>	Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Half Rails
<b>US5000-RFPKGQR</b>	Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Quarter Rails
<b>US5000-RFPKGC</b>	Semi-Electric Bed with 1633-Rolled Foam Mattress and Chrome Full Rails
<b>US5000-RPKG</b>	Semi-Electric Bed with No Mattress and FDA Full Rails
<b>US5000-RPKGHR</b>	Semi-Electric Bed with No Mattress and FDA Half Rails
<b>US5000-RPKGQR</b>	Semi-Electric Bed with No Mattress and FDA Quarter Rails
<b>US5000-SFPKG</b>	Semi-Electric Bed with 1633-Standard Foam Mattress & FDA Full Rails
<b>US5000-SFPKGHR</b>	Semi-Electric Bed with 1633-Standard Foam Mattress & FDA Half Rails
<b>US5000-SFPKGQR</b>	Semi-Electric Bed with 1633-Standard Foam Mattress & FDA Quarter Rails
<b>US5000-SFPKGC</b>	Semi-Electric Bed with 1633-Standard Foam Mattress & Chrome Full Rails