

EXTENDED CARE BEDS

# MATRIX SERIES

The Matrix bed features sealed ball bearings in the joints for smooth, quiet movement and SilverSolutions in the powder coating to provide antimicrobial surface protection.

Basic American by Graham-Field is a leading manufacturer of beds serving the senior living industry for 40+ years. We proudly manufacture our beds in the USA, in Fond du Lac, Wisconsin, which enables faster response to customers' needs.

Innovation with industry-leading technology is designed into the Basic American Zenith, Matrix, and Liberty models to provide feature-rich beds and value for every budget.



**6100**  
MATRIX BED  
*shown in chair position with  
optional two-pedal locking system*

## MATRIX SERIES

### FEATURES AND ADVANTAGES:

#### Made in the USA

- Minimal lead times
- Quality craftsmanship and manufacturing

#### SilverSolutions

- Located in the powder-coating which provides antimicrobial surface protection
- Premium Products

#### Sealed Ball Bearings

- Allows seamless, continuous, and noiseless movement
- Maintenance Free

#### Expandable Deck Available

- Pin-on wide kits extend the deck width to 39" or 42"

#### ANSI/AAMI ES60601 Safety Certification

- Tip stability testing
- Load capacity testing (twice the operating capacity)
- Roll clearance testing
- Splash proof electronics testing

#### Optional 'Plug and Play' digital electronics enable future feature additions

- Staff Control for Trendelenburg and Reverse Trendelenburg
- Chair Position
- Underbed Light
- USB Charger

#### Certified as a transport device for the CMS Emergency Preparedness Plan

- Easy maneuvering/adjustability through doorways and hallways during emergency or evacuation

#### Optional Intuitive two-pedal locking system

- Centrally located
- Easily visible and accessible at ALL heights and widths

#### Industry-leading, energy-saving electronics

- Consume up to 66% less power than transformer electronics

#### Grounded bed with right-angle hospital grade plug

- Provides added safety for both staff and residents

#### US Patented design of high/low mechanism

- Housed with the mainframe tubes; no holes are located on the top side that will harbor bacteria

#### Grid sleep surface with cross braces

- Sturdy, no-sag sleeping surface



**5100**  
MATRIX BED

USB charger (right)  
underbed light (below)



**6100**  
MATRIX BED

shown in lowest position with optional head and foot boards

### Warranty:

180	Main Frame & Welds <sup>1</sup>
48	Control Box & Actuator Motors
36	Hand Control Pendant
36	Staff Control & Cabling
12	Head Board & Foot Board
24	All Other Durable Components not otherwise listed


<sup>1</sup> 15 year frame warranty



**6100**

MATRIX BED

shown in the Trendelenburg position and with optional pedal lock

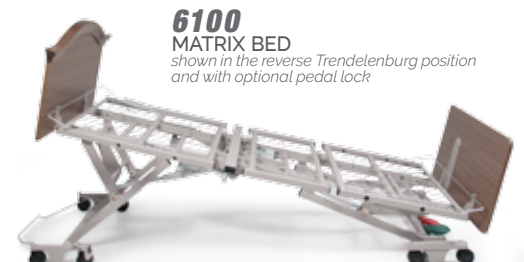

**6100**

MATRIX BED

shown in the reverse Trendelenburg position and with optional pedal lock

**WARNING: Cancer and Reproductive Harm** - www.P05Warnings.ca.gov



GF Health Products, Inc. is an  
ISO 13485:2016  
MDSAP Certified Company



+1 770.368.4700

Information contained herein is subject to change without notice.  
The most current and complete product information can be found on our website.  
[www.grahamfield.com](http://www.grahamfield.com)

