


An Inside **Look** At **LUMEX**[®] **Wheelchair Cushions**



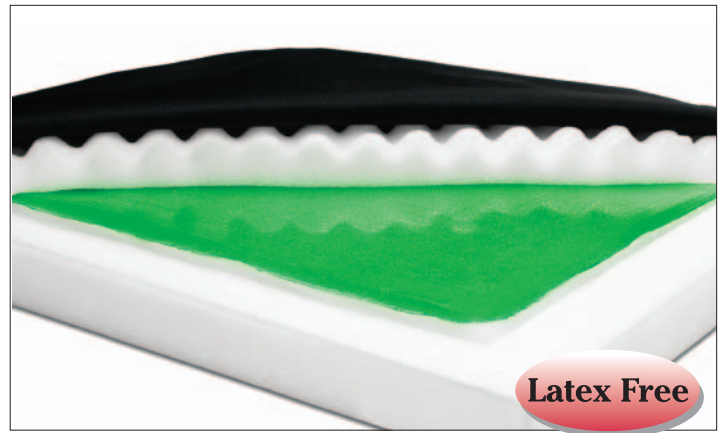

Everyday Wheelchair Cushions

HCPCS Code: E2601

For use with low risk patients

Three layers provide pressure distribution and weight equalization:

1. High density foam base

2. Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points

- Proprietary gel is a flexible solid that will not leak
- Gel is placed at three pressure zones

3. Medium density foam topper

Cover features

- Fluid-proof 600 denier cover
- Antimicrobial
- Flame retardant
- Removable and washable

2" Cushion

250 lb maximum weight capacity
18 month limited warranty

E2601

8265166	16"x 16"x 2"
8265186	18"x 16"x 2"
8265188	18"x 18"x 2"
8265206	20"x 16"x 2"
8265208	20"x 18"x 2"

General Use Cushion

HCPCS Code: E2601/E2602

For use with low-risk patients

Three layers provide pressure distribution and weight equalization:

1. High density foam base

2. Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points

- Proprietary gel is a flexible solid that will not leak
- Gel is placed at three pressure zones

3. Convoluted foam topper

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
- Antimicrobial
- Flame retardant
- Removable and washable
- Non-skid base

3" Cushion

Strap to secure cushion firmly to the chair
300 lb maximum weight capacity
2-year limited warranty


E2601

8565166	16"x 16"x 3"
8565186	18"x 16"x 3"
8565188	18"x 18"x 3"
8565206	20"x 16"x 3"
8565208	20"x 18"x 3"


E2602

8565228	22"x18"x 3"
8565248	24"x18"x 3"



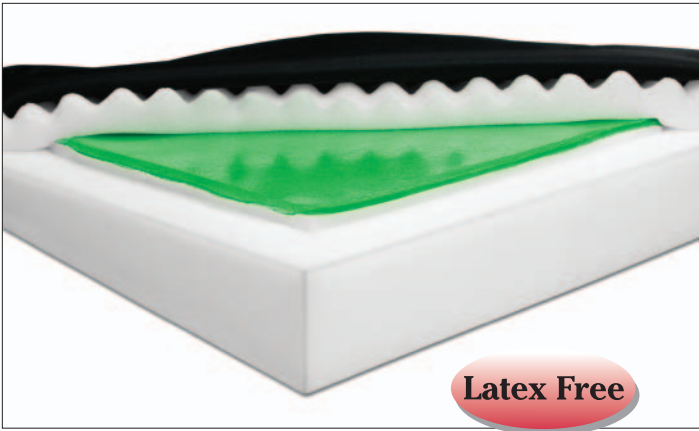

Latex Free

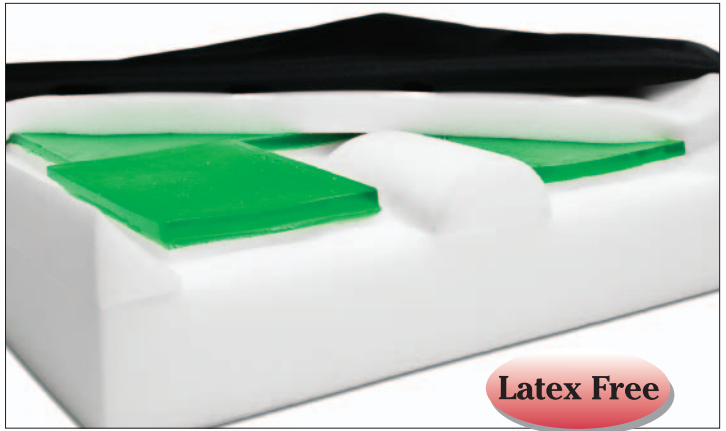

Latex Free

Basic Wedge Cushion

HCPCS Code: E2601/E2602

Three layers provide pressure distribution and weight equalization:

1. High density foam base
2. Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points
 - Proprietary gel is a flexible solid that will not leak
 - Gel is placed at three pressure zones
3. Convoluted foam topper

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
- Antimicrobial
- Flame retardant
- Removable and washable
- Non-skid base

4" in front tapers to 3 1/4" in rear
 Strap to secure cushion firmly to the chair
 350 lb maximum weight capacity
 2-year limited warranty

E2601

8220166 16"x 16"x 3 1/4 - 4"
 8220186 18"x 16"x 3 1/4 - 4"
 8220188 18"x 18"x 3 1/4 - 4"
 8220206 20"x 16"x 3 1/4 - 4"
 8220208 20"x 18"x 3 1/4 - 4"

E2602

8220228 22"x 18"x 3 1/4 - 4"
 8220248 24"x 18"x 3 1/4 - 4"

Skin Protection and Positioning Cushion

HCPCS Code: E2607/E2608

For Pressure Reduction and Positioning
 For use with moderate to high-risk patients

Three layers provide pressure distribution and weight equalization:

1. Medium density foam base
2. Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points
 - Proprietary gel is a flexible solid that will not leak
 - Gel is placed at three pressure zones
3. Visco-elastic memory foam topper for superior comfort and pressure relief

Three positioning features:

1. Side wedges provide lateral stability for proper pelvic support
2. Center pommel provides medial and lateral thigh positioning and prevents hip rotation and abrasion
3. Moderate slope from front to back so user does not fall forward

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
- Antimicrobial
- Flame retardant
- Removable and washable
- Non-skid base

4 1/2" Cushion



Strap to secure cushion firmly to the chair
 350 lb maximum weight capacity
 2-year limited warranty

E2607

8920166 16"x 16"x 4 1/2"
 8920186 18"x 16"x 4 1/2"
 8920188 18"x 18"x 4 1/2"
 8920206 20"x 16"x 4 1/2"
 8920208 20"x 18"x 4 1/2"

E2608

8920228 22"x 18"x 4 1/2"
 8920248 24"x 18"x 4 1/2"

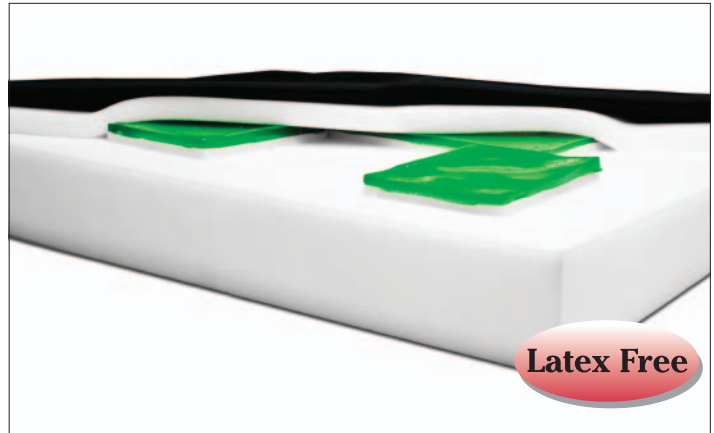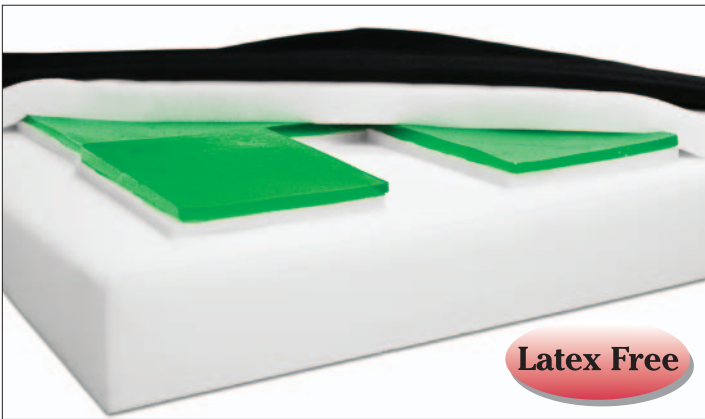

Skin Protection Cushion

HCPCS Code: E2603/E2604

For Pressure Reduction
Designed for moderate risk users

Three layers provide pressure distribution and weight equalization:

1. *Medium density foam base*
2. *Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points*
 - *Proprietary gel is a flexible solid that will not leak*
 - *Gel is placed at three pressure zones*
3. *Visco elastic memory foam topper for superior comfort and pressure relief*
Positions the user with a moderate slope from front to back

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
- Antimicrobial
- Flame retardant
- Removable and washable
- Non-skid base

3" Cushion

Strap to secure cushion firmly to the chair
300 lb maximum weight capacity
18-month limited warranty

E2603

8930166 16"x 16"x 3"
8930186 18"x 16"x 3"
8930188 18"x 18"x 3"
8930206 20"x 16"x 3"
8930208 20"x 18"x 3"

E2604

8930228 22"x 18"x 3"
8930220 22"x 20"x 3"
8930248 24"x 18"x 3"

Bariatric Skin Protection

HCPCS Code: E2604

For Pressure Reduction

Three layers provide pressure distribution and weight equalization to meet the unique demands of the bariatric patient:

1. *Medium density foam base*
2. *Non-aqueous specialized visco-elastic gel core layer reduces tissue interface pressure at critical points*
 - *Proprietary gel is a flexible solid that will not leak*
 - *Gel is placed at three pressure zones*
3. *Visco elastic memory foam topper for superior comfort and pressure relief*
Positions the user with a moderate slope from front to back

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
- Antimicrobial
- Flame retardant
- Removable and washable
- Non-skid base

4" Cushion

Strap to secure cushion firmly to the chair
650 lb maximum weight capacity
18-month limited warranty

E2604

20252218B 22"x 18"x 4"
20252418B 24"x 18"x 4"
20252619B 26"x 19"x 4"
20253019B 30"x 19"x 4"



Latex Free

Lumbar Cushions

HCPCS Code: E2611-/2612

- Contoured lumbar support provides back support and comfort
- Improves seating posture
- High density foam

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
 - Antimicrobial
 - Flame retardant
 - Removable and washable
- Back strap secures in place
 - Warranty 18 months

Designed to meet the requirements of HCPCS Code:


E2611

- 821617 16"W x 17"H
- 821817 18"W x 17"H
- 822019 20"W x 19"H

Designed to meet the requirements of HCPCS Code:

E2612

- 822219 22"W x 19"H

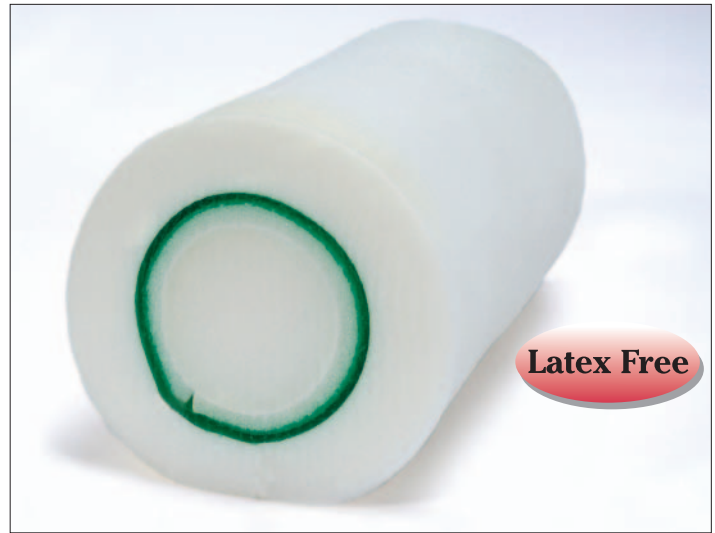

Latex Free

Lumbar Support/Neck Roll

Three layers provides support and comfort to neck or lumbar area

Jelly roll design:

1. Inner high density foam cylinder core

2. Non-aqueous specialized visco-elastic gel layer surrounding core, reduces tissue interface pressure at critical points

- Proprietary gel is a flexible solid that will not leak

3. Surrounded by 1.5" of visco-elastic memory foam

Cover features

- Fluid-proof stretch top cover, reduces shear/friction
 - Antimicrobial
 - Flame retardant
 - Removable and washable
- End handle and convenient strap provide easy handling and positioning
 - Warranty 18 months

8200 16"L x 6"D cylinder



GF Health Products, Inc. Atlanta, Georgia 30360
800.347.5678 www.grahamfield.com